PS.R2 PRECISION BALANCES









release date 05-03-2015





PS.R2 series balances represent a new standard of precision balances. They feature a new, readable LCD display which allows a clearer presentation of the weighing result. Besides, the display has a new text information line allowing to show additional messages and data, e.g. product name or tare value. New PS.R2 balances, like previously designed PS series

balances, have pans in two possible dimensions: 128x128 mm or 195x195 mm. balances with a smaller pan have a draft shield. The balance precision and the measurement accuracy is assured by automatic internal adjustment, which takes into consideration temperature changes and time flow.

PS.R2 balances feature several communication interfaces: 2 x RS 232, type A USB, type B USB and optional WiFi. The housing is made of plastic, and the pan is made of stainless steel. The balances have a possibility to weigh products out of the pan (under hook weighing) - the load hangs under the pan.

The information system is based on 5 databases, which allows for

several users to work with several products databases, and the

registered weighing results can be subject to further analysis.

DATABASES IN R SERIES BALANCES

The data is registered in 5 databases:

- weighments (up to 5000 weighments),

-ALIBI memory (up to 100 000 weighments).

- products (up to 1000 products),

databases using USB pen drives.

-users (up to 10 users),

-tares (up to 100 tares),



Parts counting





Percent



Statistics



Animal weighing



Statistical Quality Control



(GLP, Filter)





GLP procedures



Newton unit measurement





Summing



Caps lock of max indication



determination




Under-hook



ALIBI Memory





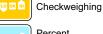

deviations







Replaceable



function



Density



weighing





Quick access to information

Direct access to functions and databases is possible from the level of keyboard.

Database - a direct access to databasis

Function - a direct access to the basic functions

F1 to F4 – programmable function and navigation keys on the menu

ALIBI memory

The used ALIBI memory is a data secure area and allows to record up to 100 000 weighment records. It ensures security of constant data register in the long time period.

There is two directions data exchange within the system thanks to

a quick USB interface. New balances allow to import and export

	PS 200/2000.R2 M	PS 210.R2 <mark>M</mark>	PS 360.R2 <mark>M</mark>	PS 600.R2 M	PS 750.R2 <mark>M</mark>	PS 1000.R2 M
Max capacity	200 / 2000 g	210 g	360 g	600 g	750 g	1000 g
Minimum load	20 mg	20 mg	20 mg	20 mg	20 mg	20 mg
Readability	1 / 10 mg	1 mg	1 mg	1 mg	1 mg	1 mg
Tare range	-2000 g	-210 g	-360 g	-600 g	-750 g	-1000 g
Repeatability *	1 / 10 mg	1 mg	1 mg	1,5 mg	1,5 mg	1,5 mg
Linearity	±2 / ±20 mg	±2 mg	±2 mg	±3 mg	±3 mg	±3 mg
Pan size	128×128 mm					
Working temperature	+10 ÷ +40 °C					
Stabilization time	2 s / 1,5 s 2 s					
Sensitivity drift	2 ppm/°C in temperature +10 ÷ +40 °C					
Interface	2 × RS 232, USB-A, USB-B, WiFi - option					
Power supply**	12 ÷ 16 V DC / 2,1 A					
Adjustment/calibration	internal (automatic)					
Display	LCD (backlit)					
Net weight/Gross weight	3,5 / 5,5 kg	3,5 / 5,5 kg	3,5 / 5,5 kg	3,5 / 5,5 kg	3,5 / 5,5 kg	3,5 / 5,5 kg
Packaging size	470×380×336 mm					


^{*} Repeatability is expressed as a standard deviation from 10 weighing cycles.

^{** 250} mA for balances without WiFi module, 350 mA for balances with installed WiFi module

Technical data:						
	PS 1200.R2	PS 2100.R2	PS 3500.R2	PS 4500.R2	PS 6000.R2	PS 6001.R2
	M	M	M	M	M	M
Max capacity	1200 g	2100 g	3500 g	4500 g	6000 g	6000 g
Minimum load	500 mg	500 mg	500 mg	500 mg	500 mg	500 mg
Readability	10 mg	10 mg	10 mg	10 mg	10 mg	100 mg
Tare range	-1200 g	-2100 g	-3500 g	-4500 g	-6000 g	-6000 g
Repeatability *	10 mg	10 mg	10 mg	10 mg	15 mg	100 mg
Linearity	±20 mg	±20 mg	±20 mg	±20 mg	±30 mg	±100 mg
Pan size	195×195 mm					
Working temperature	+10 ÷ +40 °C					
Stabilization time	1,5 s					
Sensitivity drift	2 ppm/°C in temperature +10 ÷ +40 °C					
Interface	2 × RS 232, USB-A, USB-B, WiFi - option					
Power supply**	12 ÷ 16 V DC / 2,1 A					
Adjustment/calibration	internal (automatic)					
Display	LCD (backlit)					
Net weight/Gross weight	3,6 / 5,1 kg	3,6 / 5,1 kg	3,6 / 5,1 kg	3,6 / 5,1 kg	3,6 / 5,1 kg	4,8 / 6,3 kg
Packaging size	470×380×336 mm					

^{*} Repeatability is expressed as a standard deviation from 10 weighing cycles.

 $^{^{\}star\star}$ 250 mA for balances without WiFi module, 350 mA for balances with installed WiFi module

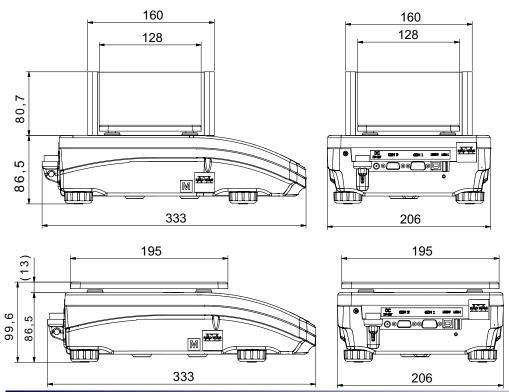

Accessories:

Antivibration table for laboratory balances	Bar code scanner USB HID			
Antivibration table SAL/STONE	Power loop output AP2-1 (plastic housing)			
Impact printer Epson	Additional LCD display "WD-6"			
Label printer Citizen	Power adapter with battery and charger ZR-02			
Printer USB PCL	Mass standard			
Density determination kit	Szafka przeciwpodmuchowa (do wag z szalką 128x128)			
Rack for under hook weighing	USB A- USB B cable (balance - computer, balance - PLC printer)			
"Tare" or "Print" foot button	Cable RS 232 (scale - computer) "P0108"			
PC keyboard USB	Cable RS 232 (scale, Epson, Citizen printer) "P0151"			
External USB memory (FAT files format)	"PW-WIN" computer software			
Bar code scanner	"RAD-KEY" computer software			

